wellbox®

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Ref. SG0906-DM-indB

www.wellbox.com

wellbox® S

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wellbox® S

You have just acquired a Wellbox®[S], medical device, for which the effects have been scientifically proven. Your Wellbox®[S] will bring you satisfaction as long as you respect the protocols mentioned in the instructional.

FIRST USE

> TREATMENT HEADS:

Open the device by pressing on the button located in front. Store the four treatment heads in the locations specifically designed for that purpose. Each head's dedicated storage area is indicated with a pictogram.

> INSTALLATION:

Insert the hose with the arrow (in relief) upwards, as indicated in the photo.

Choose one of the heads (LIFT or ROLL) and then connect it to the hose by lining up the arrows on the head and hose.

> POWER SUPPLY:

Insert the male jack into the female socket located at the rear of the Wellbox®[S]. Connect the power supply to a wall outlet by choosing first, the one adapted to your country.

> STARTING UP:

After connecting the power, press the key to turn on the device. The control panel lights up.

The Wellbox®[S] Routine, selected by the treatment menu will inform you, all the power settings to adjust.

- **The suction power (Power)** must be adjusted depending on how you feel and on your skin quality.
- **Suction sequentiality (Select)** is the number of aspirations per second.

continue suction

8 micro-aspirations per second

4 micro-aspirations per second

16 micro-aspirations per second

• The duration of the treatment (Time).

These two pictograms Tête **ROLL** and Tête **LIFT** are provided for information purposes only. They will guide you through the routine and indicate which treatment head to use according to the suction sequentiality you have selected.

• Treatments:

Select the treatment heads (LIFT or ROLL) to use, and then adjust the power settings (suction power, suction sequentiality, The duration of the treatment) according to the treatment menu.

FACE TREATMENT

For better tissue grab, move the LIFT head following a dotted line (with small hops, stopping 2-3 seconds on each point) while respecting the direction indicated by the white arrow.


The red lines indicate how to position the flaps when moving the treatment head. They must always be perpendicular to the wrinkle.

BODY TREATMENT

The white arrows indicate the direction that the treatment head must follow. When using the ROLL(3) head, use the trigger to reverse the direction of the rollers.

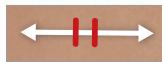The red lines indicate how to position the rollers when moving the treatment head.

Press to start the treatment. You can interrupt the treatment at any time and resume when you want by pressing .

➤ STORAGE:

To close your Wellbox®[S] after use:

- Connect the ROLL(3) head to the hose.
- Roll the hose into the tube designed for this purpose (two times around).
- Place the ROLL(3) head in the location next to the control panel.

wellbox® S

INTENDED USE

- The Wellbox®[S] allows the treatment of cellulite by smoothing the skin (reduction of the orange peel effect) visible action, according to 20 MHz ultrasound assessment, by a significant reduction in the length of the interface between the dermis and the hypodermis. Furthermore, the perimeter of the treated area is decreased, for example, a reduction of up to 3 cm is possible on a thigh which has been treated.
- The microcirculation of treated patients is improved by up to 200%, as the results recorded on the lower limbs by the laser Doppler prove it, which can be associated with a feeling of lighter legs.
- Pain such as Delayed Onset Muscle Soreness and muscular tension is reduced; following an intense effort, a 75% recovery of muscular strength deficiency is gained from the second day of treatment. Muscular aches felt by patients suffering from fibromyalgia have been reduced by 50% after 15 sessions of treatment, proven by a specific study.
- The application of Wellbox®[S] also allows the mechanical softening of the skin. All these results can vary depending on your health, age, medical history and life style.

In order to optimize the effectiveness of the Wellbox®[S], read the installation and user guide. If you have a doubt, do not hesitate to seek medical advice from your physician.

Identification label

Description of symbols

	This icon indicates that "Your device contains a part that comes in contact with the patient that is electrically isolated from all other parts of the device."
	Wellbox®[S] is marked CE 0088 as a medical device under Appendix II of the 93/42/EC directive (consolidated version including the 2007/47/EC directive).
	Carefully read the accompanying document before using the device.

› ELECTROMAGNETIC COMPATIBILITY:

TABLE 1: DIRECTIVES AND MANUFACTURER DECLARATION – ELECTROMAGNETIC EMISSIONS		
Wellbox®[S] is intended for use in the electromagnetic environment specified below. The Wellbox customer or user should ensure that it is used in such an environment.		
Emissions test	Conformity	Electromagnetic Environment - Directives
RF emissions CISPR 11	Group 1	Wellbox®[S] uses RF energy only for its internal functions. Therefore, its RF emissions are very low and unlikely to cause interference in nearby electronic devices.
RF emissions CISPR 11	Class B	Wellbox®[S] may be used in all establishments, including domestic sites and sites that are directly connected to the low voltage public power grid, which supplies domestic buildings.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations and flicker IEC 61000-3-3	Conforms	

TABLE 2: DIRECTIVES AND MANUFACTURER DECLARATION – ELECTROMAGNETIC IMMUNITY – EMISSIONS TESTS			
Wellbox®[S] is intended for use in the electromagnetic environment specified below. The Wellbox customer or user should ensure that it is used in such an environment.			
Emissions test	Conformity	Conformity	Electromagnetic Environment - Directives
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV on contact ± 8 kV in the air	± 6 kV on contact ± 8 kV in the air	Flooring should be made of wood, concrete, or ceramic tile. If the floor is covered with synthetic material, relative humidity should be at least 30%.
Fast transient/burst IEC 61000-4-4	± 2 kV for electrical power lines ± 1 kV for input/output lines	± 2 kV for electrical power lines ± 1 kV for input/output lines	The quality of the electrical power network should be equivalent to that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV between phases ± 2 kV between phase and earth	± 1 kV between phases ± 2 kV between phase and earth	The quality of the electrical power network should be equivalent to that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations IEC 61000-4-11	<5% UT [→95% voltage dip] during 0.5 cycle 40% UT [60% UT voltage dip] during 5 cycles 70% UT [30 % UT voltage dip] during 25 cycles ←5% UT [→95% voltage dip] during 5s	<5% UT [→95% voltage dip] during 0.5 cycle 40 % UT [60% UT voltage dip] during 5 cycles 70% UT [30 % UT voltage dip] during 25 cycles ←5% UT [→95% voltage dip] during 5s	The quality of the electrical power network should be equivalent to that of a typical commercial or hospital environment. If the Wellbox user requires operation to continue during power outages, the Wellbox should be powered through an uninterruptible power supply or a battery.
Power frequency magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should have levels that are characteristic of a typical hospital or commercial environment.

NOTE: UT is the AC voltage prior to the application of the test level.

› ELECTROMAGNETIC COMPATIBILITY:

TABLE 3: DIRECTIVES AND MANUFACTURER DECLARATION – ELECTROMAGNETIC IMMUNITY – IMMUNITY TESTS

Wellbox®[S] is intended for use in the electromagnetic environment specified below. The Wellbox®[S] customer or user should ensure that it is used in such an environment.

Immunity Test	Test level according to IEC 60601	Conformity level	Electromagnetic environment - directives
Conducted RF disturbances IEC 61000-4-6	3 Veff 150 kHz to 80 MHz	3V	Portable and mobile RF communications devices should not be used closer to any part of the Wellbox®[S], including its cables, than the recommended separation distance, calculated from the applicable equation on transmitter frequency. Recommended separation distance: $d = 1.2\sqrt{P}$ $d = 1.2\sqrt{P} \text{ 80 MHz at 800MHz}$ $d = 2.3\sqrt{P} \text{ 800 MHz at 2.5 GHz}$ where P is the transmitter's maximum power output in watts [W], according to the transmitter's manufacturer, and d is the recommended separation distance in meters [m]. The field strength for fixed RF transmitters, as determined by an electromagnetic survey at site ^a , should be less than the conformity level for each range of frequencies ^b : Interference may be produced near the device marked with the following symbol:
Perturbations RF rayonnées CEI 61000-4-3	3 V/m De 80 MHz to 2,5 GHz	3 V/m	

NOTE 1: At 80 MHz and 800 MHz, the highest frequency range applies. NOTE 2: These directives may not apply in all situations.
Electromagnetic propagation is affected by absorption and by reflections of structures, objects, and people.

^a The field strength of fixed transmitters, such as radio/telephone (cellular/wireless) base stations and land mobile radios, amateur radios, AM and FM radio broadcasting, and TV broadcasting cannot be theoretically predicted with accuracy. To measure the electromagnetic environment due to fixed RF transmitters, consider performing an electromagnetic survey of the site. If the field strength, measured where the Wellbox is used, exceeds the applicable RF conformity level above, watch the Wellbox to ensure that it is operating as normal. If there is any abnormal performance, additional measures may be necessary, such as reorienting or repositioning the Wellbox.

^b For the frequency range of 150 kHz to 80 MHz, field strength should be less than 3 V/m.

TABLE 4: RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATION DEVICES AND WELLBOX®[S]

Wellbox®[S] is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The Wellbox®[S] user or customer can prevent electromagnetic interference by maintaining a minimal distance between the portable or mobile RF communications device (transmitter) and the Wellbox®[S], as recommended below, based on the maximum transmission power of the communication device.

Maximum rated output of the transmitter W	Separation distance according to the transmitter frequency m		
	d=1,2VP	d=1,2VP	d = 2,3VP
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters whose maximum rated output is not shown above, the recommended separation distance d in meters [m] can be estimated using the applicable equation for the transmitter frequency, where P is the transmitter's maximum output in watts [W], according to its manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the highest frequency range applies.

NOTE 2: These directives may not apply in all situations.

Electromagnetic propagation is affected by absorption and by reflections of structures, objects, and people.

NOTES . . .

wellbox® S

WELLBOX® [S] ROUTINES

- For the best use of your Wellbox®[S] and to meet all your beauty objectives, treatment protocols (for face and body) have been specially developed by LPG® Training service.
You will find all your treatments in your private account at www.wellbox.com and on video on the «Wellbox by LPG®» Youtube page:

<https://www.youtube.com/channel/UC7gKtLhOKS308kzvrYQwvXQ>

USER GUIDE PRODUCT WARRANTY

- To reduce the environmental impact while still maintaining the level of safety and ease of use, the complete user guide is available in electronic format at www.wellbox.com.
- This guide must be read before you use your device. If you have the slightest doubt about the operation or use of your device, or if you wish to receive these instructions on paper, do not hesitate to contact

Wellbox® Customer Service hotline


00 800 699 355 269 (*toll-free in Europe only*)

or +33(0)4 75 78 68 99.

Also, take the time
to activate your warranty
directly on line:

<http://warranty.wellbox.com>

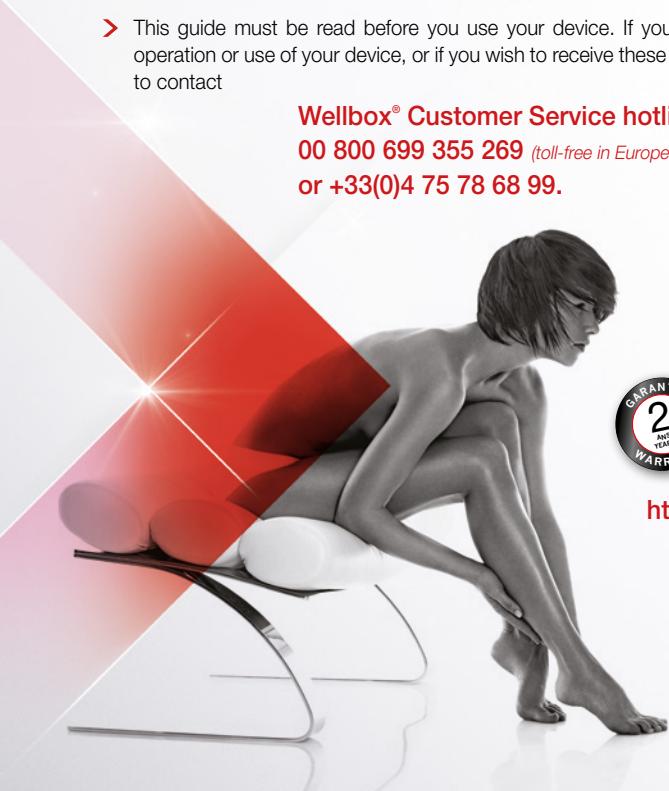CONTRA-INDICATIONS

Wellbox®[S] can be used on all healthy external parts of the body, with the exception of the eyes, ears, glottis, and genitals. The device should not be used on mucous membranes, nipples, the spine or on any suspicious moles without medical advice.

Do not use Wellbox®[S] in cases of:

- Epilepsy problems
- Blood clotting disorders
- Heart problems
- Pain or red blotches
- Tingling in the limbs
- Pregnancy
- Phlebitis (blood clots), varicose veins, or if you have undergone an operation on the veins within the last 2 months
- Fever, infection, inflammation, acute dermatitis, open wounds, recent scars, or any serious and/or progressive illness
- If taking long-term medications (anti-inflammatory drugs, anticoagulants, antibiotics or antihistamines)
- Wellbox cannot be applied on areas that are swollen, or affected by rosacea

Stop the treatment in the event of unexplained calve pain. If any side effects linked to Wellbox®[S] treatments appear during or after use, consult your doctor and contact LPG® Customer Service Department.

If you have the slightest doubt concerning the compatibility of your state of health with the use of Wellbox®[S], seek your doctor's advice.


wellbox® S

IMPORTANT SAFETY INSTRUCTIONS

- Check that the network voltage complies with the range indicated on the power supply of the device (100-240 VAC). It should be connected to a wall socket that complies with the electrical standards in force. Avoid connecting Wellbox®[S] to a multi-socket power bar. After every use and before cleaning, disconnect the Wellbox®[S] power system.
- Do not use Wellbox®[S] while taking a bath or shower. Do not place Wellbox®[S] on a wet or unstable surface, or near a bathtub, shower or sink. Use Wellbox®[S] in a dry area. Avoid splashing any water on your device. Do not touch a device that may have fallen into water. Immediately disconnect from the power supply. Do not put the hose in water. Do not use Wellbox®[S] outdoors.
- Do not drop the device, or insert objects in any openings. Do not suck up any solid bodies or liquids with Wellbox®[S].

WARNING: In order to reduce any risk of burns, fire, electric shock or injuries:

- Do not use the device under any kind of blanket or cover. Excessive heat may be produced which may provoke fire or an electric shock, or injure someone.
- Never use the device if a cable or connection is damaged, if it doesn't work properly, if it has been dropped or damaged, or if it has fallen into water. In these cases, please contact LPG® Customer Service.
- Do not pick up or move the device by its cover, cables or hose, or use them as carrying handles. Do not drop Wellbox®[S]. Do not shake or hit the device.
- Do not use Wellbox®[S] without first connecting the hose to the heads. Keep the cables and the hose away from any source of heat. Do not use the device when the temperature is above 40°C.
- Do not use the device when aerosol sprays have been used or when oxygen has been administered.
- The device must not be left unsupervised when connected. To turn the device off, press the on/off push-button, then disconnect the AC/DC power supply from the wall socket.
- Do not open or dismantle the device.

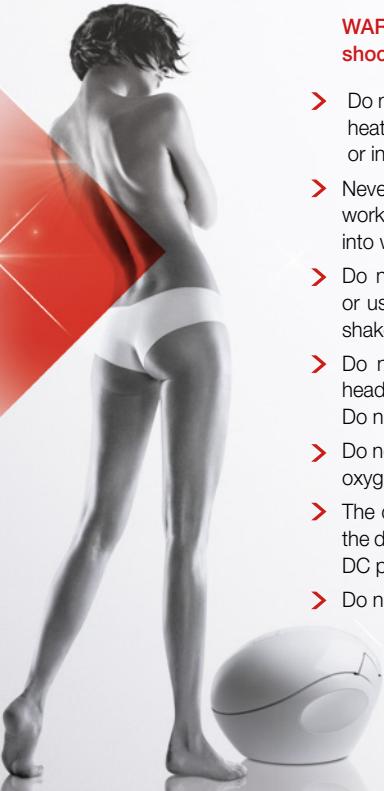SAFETY

PRECAUTIONS

- Do not use Wellbox®[S] on a child under the age of 7, and do not let a child under the age of 11 use Wellbox®[S] alone. Remember to reduce the intensity when treating a child.
- Close supervision is necessary when the device is used by, on or near children, or disabled or handicapped persons.
- Using Wellbox®[S] should never be painful or leave marks; any redness should disappear after a few seconds and there should be no marks.
- Never use Wellbox®[S] without its protection filter.
- Use the device only as described in this guide. Use only the treatment heads and accessories supplied by LPG®.
- Do not use detergents, solvents, chemical products or abrasive sponges to clean Wellbox®[S]. Use the cleaning wipes supplied, or wipes soaked in a bactericidal and fungicidal solution.
- Never use Wellbox®[S] on jewellery or piercings.
- Do not operate the device at an altitude above 2,500 m.
- It is recommended that the device be paused for 30 minutes between two 30-minute sessions when being used at a temperature above 25°C. Using Wellbox®[S] regularly makes you feel better; it is recommended that its use not be interrupted for too long a period. If the instructions and recommendations for use are not followed, you risk not achieving your desired results. Do not use Wellbox®[S] for longer than the recommended time.

Results can vary according to age, lifestyle (sports activity, diet, etc.) and morphology.

Caution: It is strictly forbidden to modify your device without the permission of LPG®.

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