User Manual


Please read the operating instructions carefully, particularly the safety precautions, before using the appliance.

ENGLISH

ILLUSTRATED GUIDE

- | | | |
|---------------------------|--------------------|----------------------------|
| 1. Control Panel | 7. Flash ready | 13. Energy Level Indicator |
| 2. Flash Button | 8. Lamp Life Light | 14. HR Lamp |
| 3. Flash Status Indicator | 9. Power Light | 15. SR Lamp |
| 4. Flash Emitter | 10. Power Button + | 16. Back Lamp |
| 5. Power Cord | 11. Power Button - | |
| 6. Operating Modes | 12. Start Buttin | |

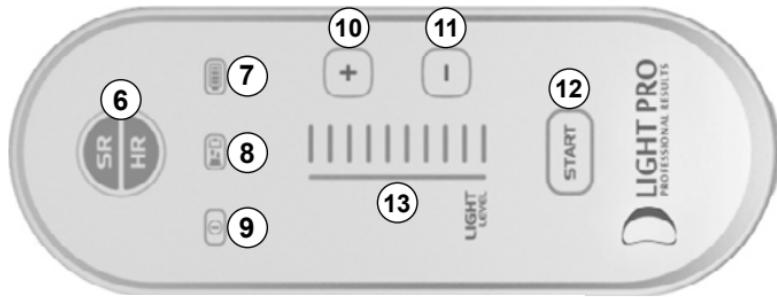Epilation and photorejuvenation - 7 cm²

LIGHT PRO
PROFESSIONAL RESULTS

20.000 spot x 2

30.000 spot x 5

USB Connection

D LIGHT PRO USER MANUAL

Dear D Light Pro Customer, thank you for purchasing this product. We are certain that you will appreciate the quality and reliability of this device, designed and manufactured by putting customer satisfaction in the foreground.

Please read the operating instructions carefully, particularly the safety precautions, before using the appliance. Keep this manual along with the relative illustrative guide, for the entire life span of the appliance, for consultation purposes. In case the appliance is transferred to a third party, also consign the entire documentation. Should certain parts of this booklet be difficult to understand or should doubts arise, contact the company before using the product, at the address indicated on the last page.

INDEX

Safety Warnings	Pag. 1-2
Contraindications	Pag. -8
Other Safety Precautions	Pag. 8-9
Installation	Pag. 10
Operating Modes	Pag. 11
How to Prepare Treatments	Pag. 12
Epilation Treatment Procedure	Pag. 14
Photorejuvenation Treatment Procedure	Pag. 15
Possible Side Effects	Pag. 16
Treatment Areas	Pag. 18
Patch Test	Pag. 19
Maintenance	Pag. 20
Problems and Solutions	Pag. 21
Technical and Environmental Specifications	Pag. 22

SAFETY WARNINGS

- This appliance is not intended for use by anyone under the age of 16; by persons with reduced physical, sensory, mental capabilities, or by inexperienced persons. Unless they have been given supervision or instruction concerning safe and proper use of the device by a person responsible for their safety. Children should not use, play with, maintain or clean the device.
- If you are using the appliance in bathroom premises, disconnect the plug from the socket after use, proximity to water is also dangerous even the appliance is switched off.
- In order to ensure better protection, it is advisable to install the power circuit feeding the bathroom a residual current device, the residual current should not exceed 30 ma (ideal 10 ma). Ask your installer for advice.

CAUTION! Do not use this appliance in the vicinity of bathtubs, showers, sinks or other containers containing water! Do not immerse the appliance water!

- After removing the appliance from the packaging, make sure that the product is complete and in perfect condition. Check the integrity of the supply, the design and any presence of transport damages. In case of doubt, contact the authorized service department.
- Before connecting the appliance, ensure that the power supply voltage corresponds with the rating quoted on the equipment identification nameplate. The nameplate is located on the appliance's power supply.
- Do not expose the appliance to moisture or atmospheric agents (rain, sun).
- Always disconnect appliance from the power supply when not in use and before cleaning or servicing.
- This appliance must be used only for the purpose for which it has been designed as a pulsed light fixture for home use. Any other use is considered improper and therefore dangerous.
- This device is not suitable for use in the presence of flammable anesthetic mixture with air or oxygen or nitrous oxide.
- Do not use the appliance with wet hands, damp or bare feet.

- Do not pull the power cord or the appliance to unplug it from the power socket.
- If the appliance fails or malfunctions, turn it off and do not tamper with it.
- The appliance can be used only on human hair. Do not use it on animals.
- Do not obstruct the ventilation holes.
- Do not use appliance if it is cracked or broken.
- Do not use the appliance if the light is damaged, or the surface treatment is cracked, chipped or missing.
- Never look directly at the light coming from the treatment area.
- Never point the treatment surface by trying to emit a light impulse from an open space.
- In case of side effects or complications, it is advisable to contact your doctor or dermatologist.
- Protective eyewear is highly recommended when using the appliance.

CONTRAINDICATIONS

Important Safety Information Read before use

This appliance is not intended to be used by anyone. Read and take into account the information included on the following paragraph before use. For more information and personal advice you can visit the website www.dlight-ipl.com

If you are unsure if the appliance is safe for you to use, please ask your doctor or dermatologist

DO NOT USE DEVICE IF ANY OF THE FOLLOWING CONDITIONS IS PRESENT AT THE TIME OF TREATMENT:

- Pregnant or breast-feeding
- If you were exposed to strong sunlight or an artificial tanning machine during the past 28 days.

- If you have a tattoo or permanent makeup on the area to be treated.
- If you have dark brown or black spots, such as large freckles, birth marks, or warts on the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sore in the area to be treated. Wait for the affected area to heal before using the appliance.
- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including anti-inflammatory agents (such as aspirin, ibuprofen, paracetamol, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulphonylureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acid (BHA), Retin-A®, Accutane® and/or topical retinoids.
- If you have abnormal skin conditions, such as vitiligo or albinism.
- If you have varicose veins in the treatment area.
- If you are currently or have recently been treated with Alpha-Hydroxy acids (AHA), Beta-Hydroxy acids (BHA), Retin-A, topical retinoids or azelaic acidH.
- If you have abnormal skin conditions, caused by diabetes, for example, or other systematic or metabolic diseases.
- If you have been treated with Accutane® (isotretinoin) within the past 6 months.
- If you have been taking on steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the area of treatment, unless you have consulted your physician and received preventive treatment before using the appliance.
- If you suffer from epilepsy.
- If you have an active implant, such as a pacemaker, an incontinence device, insulin pump, etc.
- If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or areas of potential skin malignancies.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated. In case of doubt, consult with your physician.

► SKIN TONE

The appliance removes unwanted hair by selectively addressing hair pigment. Varied quantities of pigment (melanin) are also present in the surrounding tissue of. The quantity of pigment in a particular person's skin, which is manifested by their complexion, determines the degree of the risk they are exposed to using any method of light-based epilation.

Before using the appliance, compare your complexion with the skin tone chart shown on the product packaging

Do not use the appliance on a naturally dark complexion!

Treating the dark skin with the appliance may cause sunburn, blisters and skin color changes (hyper or Hypo-pigmentation).

► TANNED SKIN

Tanned skin particularly after sun exposure, contains large amounts of melanin. This applies to all skin types and all complexion, including those who do not seem to tan quickly. The presence of large quantities of melanin exposes the skin to higher risk when using any method of light based treatment.

Avoid exposure to sun and/or tanning lights for 4 weeks before using the appliance!

► EYE CAUTION

Never use the appliance around or near the eyes!
The appliance can be used for the treatment of facial area. Extra caution should be taken during the treatment to avoid eye area.

► NIPPLES GENITALS AND LIPS

Do not use the appliance on nipples or genital organs (male or female) and on the lips.

OTHER SAFETY PRECAUTIONS

ENERGY LEVEL

Choose your energy levels CAREFULLY! Energy level refers to the light pulse that is projected on the skin during use, from the lower level to the highest level. Energy level indicator light (13) indicates the energy level to be used. As the energy level increases, it improves results as well as the risk of side effects (See "Possible Side Effects" below).

Always begin treatment using the lowest energy level!

Only if you experience little or no discomfort during and after use of the device energy at the lowest energy setting, you can increase the energy level, and so on for each subsequent session.
For detailed instructions on energy level setting, see "Inspection Before Use".

- Always shave the area to be treated and make sure the skin is clean and dry before using the appliance.
- Cover tattoos, lips and nipples before applying the appliance.
- Cover dark brown or black spots, such as large freckles, birth marks, moles or warts before using the appliance.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol), or acetone to clean the skin before using the appliance.
- Do not use the appliance on any body part where you might later want hair.

► PATCH TEST

Before you begin the treatment, it is recommended that you perform a patch test to identify your skin's highest level of tolerance to the pulsed light. For instructions on how to perform a patch test, see "Inspection before use" section.

► DURING TREATMENT


Not treat the same area of skin more than once per session! Do not repeat the spots!
HR spot = 2 flash **SR spot = 5 flash**


- If blisters are formed on the skin, or if the skin burns, IMMEDIATELY STOP USING THE APPLIANCE.
- Always make sure that the treatment surface points toward the skin, and is completely in contact with it.
- Remove the surface treatment from the skin if it becomes excessively hot.

INSTALLATION

1) Connect the Power Cord (5) See illustration below.

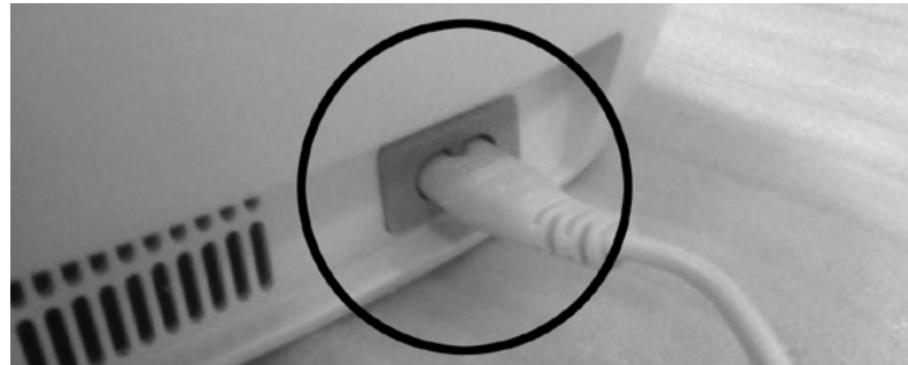2) Mounting Lamps (HR and SR)

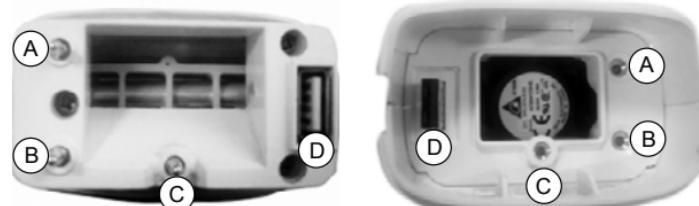Pay attention when inserting the lamp. Behind the lamp there are 3 metallic pins (a b c) and a USB socket (d), during the installation of the lamp, be careful to match these elements, equivalent to the handpiece. The installation and removal of the lamp does not require any equipment and is done by hand. N.B. The USB socket is used only to connect the lamp to the machine, any other use may be dangerous and damage the unit!

OPERATING MODES

Depending on the type of lamp mounted, the appliance automatically selects the correct mode Operation.

- HR for Epilation Mode
- SR for Photorejuvenation Mode

► HR MODE (EPILATION)

In HR mode, when the Flash Button (2) is pressed, 2 consecutive flashes are emitted with a range of 0.5 sec, the first flash(weaker) will serve to pre-heat the skin and bulb preparing to receive the second flash (more powerful) that will effectively neutralize the matrix. The two flashes that are emitted at each button push are used to treat the same area, it is recommended that you do not move the flash window during this phase.

► SR MODE (PHOTOREJUVENATION)

In SR mode, when the Flash button is pressed (2), 5 consecutive flashes are emitted with a range of 0.5 sec. In such a way to equally divide the power of the spot, and stimulate the production of collagen in a more profound and prolonged manner, without excessively stressing the skin.

► OPERATING THE HANDPIECE

There are lights (3) on the handpiece which indicate the status "Ready for Flash". When the lights are green and blink in sequence, it means that the flash is loading, when the lights become fixed blue, it means that D light Pro is ready for flash.

► HOW TO PREPARE THE TREATMENTS

Before starting treatments, prepare the treatment area by thoroughly cleaning the skin and remove makeup. Hair must be shaved with a razor or a hair clipper, it is preferable to use a clipper as you immediately obtain an optimal length of 1mm, and maximizes the absorption of pulsed light (in epilation treatment), using the clipper will also avoid any skin irritation.

Given that the pulsed light is absorbed by dark colors, it is necessary to protect small, dark birthmarks and tattoos, cover it with white eyepencil makeup, or with a specific product like the D Light Patch Pro.

Do not let the pulsed light pass on large or irregular areas! In the presence of doubt, ask your dermatologist before initiating treatments!

Do not use the appliance on naturally dark complexion or on tanned skin!

Do not use the pulsed light on varicose veins and any type of skin lesions.

► SELECTING THE ENERGY OF TREATMENT

It is always advisable to start from a low energy level and gradually increase to find the right intensity according to your own tolerance. The energy level is selected by pressing the + - buttons (10 and 11), the selected energy level is on the level indicator (13).

The following are intensity levels recommended for epilation and photorejuvenation:

- Energy from 4 to 7 for lighter skins.
- Energy from 3 to 5 for darker skins.

As the energy level increases, so does the results, as well as the risk of possible side effects, you can push to power 10 gradually and only if you do not feel discomfort.

If blistering occurs on the skin, or if you feel burning, IMMEDIATELY STOP using the device.

► EPILATION TREATMENT

During epilation treatment, it is important to follow a pattern that allows you to cover the entire area to be treated, avoiding double pass or overexposures. Once the area to be treated is determined, point the Flash window (4) in contact with the skin and press the Flash button (2). When the Flash button is pressed, the LEDs on the Handpiece (3) become green, then 2 consecutive flashes that serve to treat the same area. Once the flashes are emitted, the LEDs on the Handpiece (3) start flashing sequentially to indicate that the new spot is being loaded. When the LEDs on the Handpiece stop flashing and become all blue, it means that the flash is loaded and ready to handle the next area.

► HAIR REMOVAL TREATMENT CYCLE

The number of treatments required to maintain a hairless skin, also depends on the individual's hair regrowth rhythm , the recommended treatment cycle are as follows:

- 0–3 months- 1 treatment every 2 weeks
- 3-6 months - 1 treatment once a month
- 6-9 months - plan your treatments according to your needs
- Over 9 months - Only 1 or 2 treatments required per year.

At the end of the epilation treatment, it is highly advisable to apply a softening cream to balance the ph such as a D Light post treatment cream.

► PHOTOREJUVENATION TREATMENT

During the photorejuvenation treatment, it is important to follow a pattern that allows you to cover the entire area to be treated, avoiding double pass or overexposures. Once the area to be treated is determined, point the Flash window (4) in contact with the skin and press the Flash button (2). When the Flash button is pressed, the LEDs on the Handpiece (3) become green, then 2 consecutive flashes that serve to treat the same area. Once the flashes are emitted, the LEDs on the Handpiece (3) start flashing sequentially to indicate that the new spot is being loaded. When the LEDs on the Handpiece stop flashing and become all blue, it means that the flash is loaded and ready to handle the next area.

► PHOTOREJUVENATION TREATMENT CYCLE

Protocol for photorejuvenation of face-neck-collar and back of the hand.

In the photorejuvenation cycle, it is advisable to schedule a session every week and continue until the desired results are achieved. If pulsed light epilation treatments are also performed on the same area, it is advisable to allow at least one week between epilation treatment and treatment of photorejuvenation.

At the end of the photorejuvenation treatment, you can apply an emollient cream and Anti-Aging Moisturizer.

AFTER TREATMENT

Do not expose the skin to the sun and/or tanning lamps after treatment. For two succeeding weeks after each treatment session, in case of sun exposure, protect the treated skin with a minimum SPF30 cream.

POSSIBLE SIDE EFFECTS

When used in accordance with the instructions, side effects and complications associated with the use of the appliance are rare. However, any cosmetic procedure, including those intended for household use, involve some risk. Consequently, it is important to understand and accept the risks and complications that may occur with designed pulsed light epilation systems for household use.

NEGATIVE OCCURRENCE	PROBABILITY OF OCCURRENCE	NEGATIVE EFFECT	PROBABILITY OF EFFECT
Multiple pulses accumulated one over the other or overlapped on the same skin area.	Minimum	Minimal Skin Discomfort	Minimum
		Skin Redness	Minimum
		Increased Skin Sensitivity	Minimum
		Wounds and Burn	Rare
		Scars	Rare
		Pigmentation Disorder	Rare
		Excessive swelling and redness	Rare
		Infection	Negligible
		Bruises	Negligible

Minimal Skin Discomfort

Despite the domestic pulsed light treatments are generally tolerated very well, most users try a annoyance during use, generally described as a slight warm sensation or a pinch on the treated skin areas. This sensation can be an alert during application or for a few minutes after application. Any sensation beyond this minimal annoyance is not normal and it means that you need to interrupt the use because your skin is not able to tolerate epilation treatment, or the energy level is too high.

Skin Redness

Skin may be reddened immediately after using appliance or during the following 24 hours. Redness generally disappears within 24 hours. Consult your physician if the redness does not disappear within 2 or 3 days.

Increased Skin Sensitivity

The skin of the treated area becomes more sensitive and can be dry or scaly.

Skin wounds and burns

Burns or wounds may occur after application. They may take a few weeks to heal, and, in very rare cases, may leave a visible permanent scar.

Scars

Although it is very rare, permanent scars may be formed. Generally when scars are formed, they have a flat lesion and white skin (hypertrophic). However, the scar may be very large and red (hypertrophic) or large and extends beyond the margins of the wound itself (keloid). Subsequent aesthetic treatments can be necessary to improve the appearance of the scar.

Pigmentation Changes

In some rare cases there is a risk of temporary pigmentation (greater pigment or Brown bleaching) or IPO-pigmentation (whitening) on the surrounding skin area leeping. The risk of skin pigmentation change is higher for people with dark skin complexion

Redness and excessive swelling

In some rare cases the skin can become very red and swollen. This is more common in sensitive areas of the body. Redness and swelling should decrease between 2 and 7 days later, it must be cured with frequent applications of ice. You can gently wipe the area, but avoid exposing to the sun.

Infection

Infection of the skin is very rare but it is a possible risk due to scalding or injury caused by the device.

Bruises

Very seldom use of the device may cause blue-violet bruises that may last from 5 to 10 days.

► EPILATION TREATMENT AREA

The appliance can be used to remove unwanted hair from the body on the following areas: Armpits, bikini line, arms, legs, back, chest. It can also be used on the feminine face

In the area under the cheekbones and without getting too close to the eyes and wearing the protective goggles provided. It is important to only treat areas of the face where superfluous hair is present visible, thick and dark. Thin blond hair on the chin area should not be treated and the area under the chin.

► PHOTOREJUVENATION TREATMENT AREA

The appliance can be used to remove unwanted hair from the body on the following areas: Armpits, bikini line, arms, legs, back, chest. It can also be used on the feminine face

In the area under the cheekbones and without getting too close to the eyes and wearing the protective goggles provided. It is important to only treat areas of the face where superfluous hair is present visible, thick and dark. Thin blond hair on the chin area should not be treated and the area under the chin.

NOTE:

- Treating the same body part more frequently than recommended does not improve results. Due to hormonal and physiological changes or other type, hair follicles, latents can become active. Hair removal sessions or occasional maintenance treatment can therefore be necessary.
- Individual reaction depends on hair type and biological factors that can affect the hair growth cycle. Compared to the recommended epilation program, some people may react faster or slower.
- The appliance is not effective on very clear hair, on naturally white, gray, light blond and red hair (see color scale on the packaging). If the hairs of your body are one of these colors, the device will not be effective.

► INSPECTION BEFORE USE**ENERGY LEVEL**

The appliance can be set to 10 energy levels, from 1 (the lowest) to 10 (the highest). Represented by 10 energy level indicator lights (13).

As the energy level increases, so does the results, as well as the risk of possible side effects. It is advisable to perform a test patch to identify the right energy level according to proper tolerance.

► HOW TO PERFORM A PATCH TEST

Because the skin tone varies from area to another, you should perform a patch test to each body area to be treated. First check your skin tone comparing it with the skin color chart on the instruction manual and packaging.

- Identify a small area to be treated.
- Shave the area to be treated.
- After having set to the lowest level, apply the handpiece to the area and emit the first pulse. If you do not feel any kind of discomfort, increase the energy to the next level, move the handpiece to the adjacent part of skin and emit a new flash.

Increase the energy level on every part of the next skin until reaching the highest level (but it must remain bearable depending on skin type). During treatment it is possible (and normal) to feel a warm sensation on the skin. If there is a feeling of discomfort, it means that the highest energy level endured the previous one.

- After 24 hours check the area: the persistence of redness indicates that the skin tolerated the highest energy level on the previous one. Although it is extremely rare, if hypersensitivity persists for more than 24 hours, it is advisable to consult a physician.

CAUTION: The energy directly delivered affects the result, but because of the principle treatment, you should not exceed the maximum threshold that your skin is able to tolerate. This threshold is normally high for light skins and low for dark skins, because the presence of melanin darkens the skin, increases heat production. Do not exceed the highest threshold that your skin can tolerate without any kind of discomfort.

MAINTANANCE

CLEANING AND CONSERVATION

It is recommended to clean the Flash window after each epilation session:

- Switch off the appliance by holding the ignition button down (9) and disconnect the Power plug (6)
- Gently wipe the flash window of the device (11) with a clean dry cloth. For proper cleaning of the flash area, it is highly recommended to use Kera Flux detergent at luce-pulsata.net
- Store the appliance in a cool dry place, away from heat and water sources
- When the light is empty, the counter reaches zero.

NOTE:

If large stains appear inside the treatment surface (11) or if the latter is cracked or broken, do not use the appliance and contact the D Light Service center. After replacing the lamp, always set the energy level to lower than the one used normally. In the next session, increase the energy level again only if the there is no abnormal discomfort during or after treatment.

Never immerse the appliance or any of its parts in water!

PROBLEMS AND SOLUTIONS

PROBLEM	CONTROL / SOLUTION
The appliance does not start	Connect to a better power socket
	Contact our customer support
The contact sensor does not detect the skin	Replace the lamp
	Contact our customer support
The appliance does not emit pulses	Check if the display does not blink when in contact with the skin
	Replace the lamp
	Contact our customer support
The device buttons are not pressure sensitive	Contact our customer support
The machine emits 3 consecutive beeps and does not shoot the flash	The problem can occur when you switch to the next treated area without touching off the skin, but slide it on the same area. To unlock the machine, point the flash window toward the vacuum and press the firing button. In this way, the contact sensor is reset and the machine resumes to operate normally.

TECHNICAL SPECIFICATIONS

Type of Lamps	HR (Epilation) SR (Photorejuvenation)
Wavelength	530nm~1100nm / 650nm~1100nm
Area spot	4.4cm*1.6cm
Energy density	0 ~ 7J/cm ²
Pulse number	2 (HR) / 5 (SR)
Energy level	1 ~ 10
Power supply	AC100V-240V,50Hz/60Hz
Rated input power	150W
Net weight	1.75kg
Equipment size	L270mm*W230mm*H120mm
Package size	L300mm*W300mm*H135mm

► ENVIRONMENTAL SPECIFICATIONS

- Operating Ambient Temperature : 10°C ~ 30°C
- Storage Ambient Temperature -10°C ~ 55°C
- Degree of Operation Humidity: ≤70%
- Degree of Moisture Storage: ≤93%
- Maximum operating altitude: 860hPa ~ 1060hPa
- Do not expose device directly to the sun or other radiation

If the appliance has been stored in a room with temperature too cold for days without being used, you will need to bring it back to 15 ° C or more before you can use.

WARRANTY CERTIFICATE

This warranty is without prejudice to any rights established for the consumer by the consumption Code (Legislative Decree No. 206/2005) and art.1519 bis and following the Civil Code. SMA S.R.L. With current headquarters in Via Bianche 336010 Carrè (vi) guarantees the buyer of the product (buyer) that under normal conditions of use the electronic appliances produced, assembled and distributed by SMA S.R.L. with "D light" brand (products) are free of material and workmanship defects and have been constructed in compliance with current regulations; The same products were thoroughly tested and subjected to stringent inspections by quality control. The present is the only valid guarantee; Nobody is authorized with Exclusion of the SMA s.r.l., to extend or modify the terms of the same, and even more to release other written or verbal warranties. In any case of product replacement, as well as the extension of the warranty, is excluded, as a result of the defect of conformity

► WARRANTY PERIOD AND COVERAGE AREA

The products-intended as new products-are guaranteed for a period of 2 years from the date of purchase, which is proven by presenting the certificate and the fiscal receipt; while purchases by invoice and URP.IVA are guaranteed for a period of 1 year from the date of purchase

pursuant to law. You do not accept purchase documents with modifications and/or additions. Failure to exhibit valid documents and proof of purchase in the above terms will result in the warranty and the charge of the entire buyer's repair costs. SMA S.R.L. does not respond to the loss of the warranty certificate and will not be liable to issue duplicates. After 2 years from the date of purchase, the product will no longer be covered by warranty and assistance may be given by the technical assistance centers (d light) authorised by the SMA S.R.L. That will charge the buyer with the current tariffs. Any form of warranty expires in any case after three years after the product has been eliminated from the official list of SMA S.R.L. Repairs made under warranty will not give any extension or renewal of the same.

► WARRANTY SERVICES

All services covered by this warranty shall be subject to the actual defect of the product pursuant to this warranty. If the damage and/or defect is not covered by this warranty and/or the exclusion of the warranty, The buyer will be charged all costs related to verification, labor, replacement parts and transport of the product. The warranty may be asserted and will be effective only with the product is presented at the D Light Centre authorised by SMA S.R.L. or the direct seller. However, the product must be repaired at one of the light centers authorized by the SMA S.R.L. Under no circumstances shall the seller or buyer of the product be able to or by third parties other than authorized light centres, repairs or replacements of the product; In the case of infringement of this prohibition, the product will not be considered in warranty. It will be the buyer's responsibility to communicate to the D Light Center or to the direct seller of the product, any defects and vices of the same and the willingness to use the guarantee. The expenses and risks of transport to and from the aforementioned D Light centers will be charged to the buyer. The intervention of the authorized D light centres, the withdrawal of the product, its repairs, and its eventual substitution will occur within the time limits that are compatible with the service to be carried out, unless repairs and/or replacements are difficult. Any defective component replaced will be owned by the SMA S.R.L.

► EXCLUSIONS FROM WARRANTY

All parts that are defective or damaged due to transport, negligence and carelessness in use, incorrect installation or maintenance performed by the buyer or by unauthorized personnel are not covered by any warranty. Improperly or otherwise from the one for which the product was constructed, i.e., due to circumstances which, however, can be traced to product manufacture defects. In any case, glass parts and accessories are excluded from the warranty. Do not prove that it is a manufacturing defect. The technical interventions inherent in the installation and/or testing of the equipment are also excluded from the product, in addition to parts subject to normal wear.

In any case, the warranty does not cover the products with:

- Serial number is unreadable, deleted modified or removed;
- Defects due to any use not complying with the detailed instructions stated in the user manual;
- Defects due to faults or fluctuations in electrical current and/or electrical system, air conditioning, Humidity or other environmental conditions;
- Defects due to damages caused by the use of accessories, products or related components which are not part of the product covered by the warranty;
- Damages and malfunctions resulting in use.

The warranty does not cover any accessories such as cartons, carrying cases, etc.

► LIMITATION OF LIABILITY

Compensation for direct or indirect damages of any kind to persons or things because of misuse or suspension of use of the product is not applied. In particular, the SMA s.r.l. declines any liability for any damages that may arise, directly or indirectly, to persons, pets or things, due to the failure to comply with all the requirements set forth in the product's user manual and concerning, above all, the warnings regarding installation, use and maintenance of the product.

INFORMATIVE EX ART. 13 URD.LGS. 196/03

The personal data communicated by the buyer in the coupon of the warranty certificate issued by the dealer or in the form submitted to the authorized service center of the authorised representative will be processed by SMA S.R.L. with a computer aid and manual tools with the sole purpose of making the buyer use this warranty and in order to facilitate the intervention of the D Light Center authorised by SMA S.R.L., as well as for managerial purposes related to the exercise of the economic activities of the SMA S.R.L. And for statistical purposes, marketing and commercial information relating to the activity carried out by SMA S.R.L. Such data may be communicated to other group of SMA S.R.L. companies for purposes related to the organization of the warranty service, assistance and maintenance. For the possible exercise of rights provided in art. 7 of Legislative decree. 196/2003, the buyer may contact SMA S.R.L.-Via Bianche 3, Carrè 36010 (vi)- the holder of the data processing. In any case, the buyer may immediately oppose the processing of his personal data for marketing and commercial information, by writing to support@dlight-ipl.com

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WARRANTY CERTIFICATE

This appliance has been submitted by the manufacturer to stringent testing according to international standards. We are therefore confident of its high quality. In the event that malfunctions occur, we invite you to contact the authorized service centers that will perform technical intervention with the use of original spare parts.

CLIENT DATA

Name and Surname

Street and Number

CAP City Province

Country

Issued by:

SMA s.r.l.

Via Bianche 3
36010 Carrè (VI)
Italy

PIVA 03547050249

English customers support
+ 39 0445 1740014

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Soporte al cliente Español
+ 34 9373 70490

Multilanguage Support
support@dlight-ipl.com

CALENDAR OF TREATMENTS

CALENDAR OF TREATMENTS

CALENDAR OF TREATMENTS

MADE IN PRC