BILANCIA DIGITALE & ANALOGICA DIGITAL & ANALOGIC SCALE

**ISTRUZIONI PER L'USO
OPERATING INSTRUCTION
INSTRUCCIONES DE USO
ИНСТРУКЦИЯ ПО ПРИМЕНЕНИЮ**

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PRODUCT FEATURES

- Analogue and Digital (LCD) display are synchronized.
- Auto switch on / Auto switch off
- Round cockpit display
- Backlight on LCD & dashboard
- Reminder function (alarm after 24 hours)
- Low battery indicator


FIRST-TIME USAGE SETUP

1. Insert three fresh AAA batteries. Note the correct position of positive (+) and negative (-) battery terminals. (Fig#1)
2. The LCD will fully display for 5 seconds then turn off
3. Tap the center of the scale once to fully display LCD and then release to activate calibration. When "0.0" display on LCD, the calibration is done.

Note: When re-inserting battery or replacing battery, please follow above step 3 to activate Zero Calibration.

Fig 1

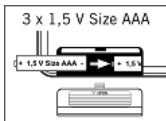WEIGHING MODE SWITCH

Adjust the weighing mode switch located on the bottom of scale to set your weighing unit (Kg/lb). (Fig#2)

Fig 2

OPERATION INSTRUCTION:

1. Step directly on the scale, and the product will automatically power on for weight measurement.
2. Keep as still as possible. The LCD display will flash 3 times and show your weight.
3. Indicated weight keeps displaying about 15 seconds when you stand on the scale. If step off the scale, weight on the display will disappear after 5 seconds to conserve the battery power.

REMINDER FUNCTION SETTING (ALARM AFTER 24 HOURS)

1. Turn-on setting: Press the button (24H/Alarm) (Fig #3), an alarm icon  will display on LCD. It will beep at the same time the next day reminding you to measure weight.
2. Turn-off setting: Press the (24H/Alarm) button again, and the alarm icon  will disappear on LCD.

Fig 3

Stop beep sound: When it beeps, just step on the scale to stop the beep sound. If you do not step on scale when beeping, the beep sound will automatically shut off after about 20 seconds.

NOTE:

1. Glass surface will be slippery when wet. Please keep surface dry and only step onto the scale when dry.
2. Never use harsh detergents or bacterial cleaners. Periodically wipe the scale surface with a DRY cloth to remove dust.
3. Always weigh yourself on the same scale placed on a hard flat surface. A hard surface is best for maximum accuracy.
4. Remove batteries if you do not use the scale for a long time.

TROUBLE SHOOTING:

In case of	Display	Action
Overloading	-----	Step off the scale
Low Battery	Lo	Replace batteries
Indicator not returning to Zero		Reset Product: Remove the battery and then insert.

PRODUCT SPECIFICATIONS:

- Loading Capacity: 180kg(396 lb)
- Graduation d=100g/0.2lb
- Workable body weight: 5 to 180kgs/ 11 to 396lbs
- Battery: AAA x 3pcs (4.5V)
- Indicator deviation: +/- 1kg (about 1.5 degree)
- Reminder function: beep alarm for 20 seconds and auto off
- LCD size: 80 x 32.5mm

CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 2 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc..
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@maniquick.com.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases:

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;

- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;
- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@maniquick.com
www.maniquick.com

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

Send this coupon in case of repairs:

Product: MQ917

Serial Number:

if available

Date of purchase:

			/				/				
--	--	--	---	--	--	--	---	--	--	--	--

Dealer's Stamp:

--

Buyer's Full Name:

Street/Square:

N°:

City and State:

Postal Code:

Country:

Phone Number:

E-mail:

@

Problem description:

.....

.....

.....

Signature:

Date:

			/				/				
--	--	--	---	--	--	--	---	--	--	--	--

WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

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DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

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SECURITY OF ALKALINE BATTERIES INTO PRODUCTS

1. Fire - Never delete the batteries throwing them in a fire. they could break and release the substances contained inside them.
2. If ingested or inserted, seek immediate medical assistance.
3. Do not store or leave the batteries in places accessible by children.
4. Remove old or weak batteries from appliances.
5. Do not mix old and new batteries in the same apparatus. There may be losses. The replacement of all batteries in a device must be done simultaneously.
6. Never mix different types of batteries in a single device. There may be losses.
7. Install carefully the batteries in the correct polarity (+ and -) according to the instructions provided by the manufacturer of the product
8. Do not leave the Products turn on with batteries when they are depleted. There may be losses.
9. Remove the batteries from any device that must remain unused for long periods of time
10. Avoid contact of the batteries with metal objects or other batteries because they may experience losses or they could break.
11. Do not recharge batteries because there may be some losses or they could break.
12. Not groped to disassemble the batteries
13. Store batteries in a cool and dry. Avoid sudden changes in temperature. Store the batteries in their original packaging before use.
14. RESPECT THE ENVIRONMENT - DISPOSE OF EXHAUSTED BATTERIES IN THE CONTAINERS

The battery and the possible Electronic Equipment -AEE- or anything that contains it must not be disposed as urban waste, but in compliance with Directive 2002/96/EC, 2006/66/EC and of National Laws transposition for the disposal of the product end of life, and must be subjected to a separate collection to avoid damage to the environment. Failure to comply with the above is sanctioned by Law.

IMPORTATO DA:
SANICO Srl
Via G.Ferraris, 37/41
20090 Cusago (MI) ITALY
Tel (+39) 02.90.39.00.38 - Fax (+39) 02.90.39.02.79
info@sanicare.it - www.sanicare.it