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LOW-FREQUENCY MAGNETOTHERAPY MAT

Globus presents a new device for magnetotherapy which can be combined with all Magnum devices.

The brand-new **MAT 100** is a mat containing four internal solenoids and enabling the treatment of an extremely large surface without placing the classic diffusers on the body.


USER-FRIENDLY: **Mat 100** can be used effortlessly everywhere: on the couch, on the bed or on a chair, by simply placing it on the body.

COMFORTABLE: Mat 100 is made up of many layers of soft material, to level the surface and for you not to feel the thickness of the internal solenoids.

MAT 100 can be used singularly or can be easily connected to other mats (**E-MAT and N-MAT**), to form a true 70-cm-wide and 176-cm-long mattress.


Instructions for use

Globus magnetotherapy mats contain four solenoids, north-oriented to the front and south-oriented to the back.

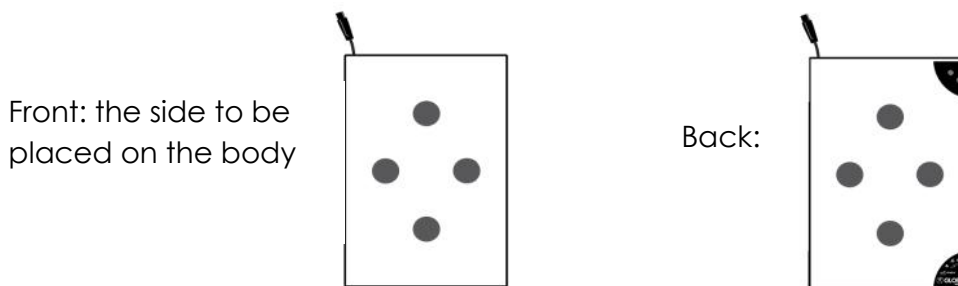
MAT 100

Mat 100 contains four internal solenoids and can be connected directly to your Magnum device through the connector.

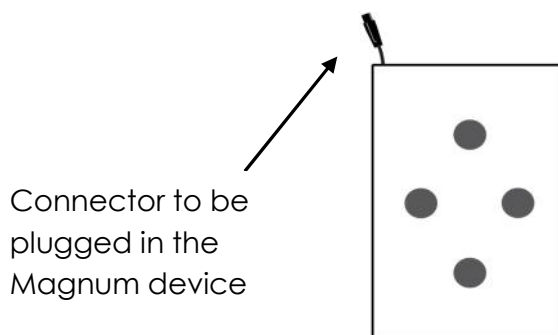
To connect Mat 100 to the device, plug the connector in the slot located in the upper part of the unit (see the picture below). When plugging in the cable, the grooves of the cable have to be oriented downward.

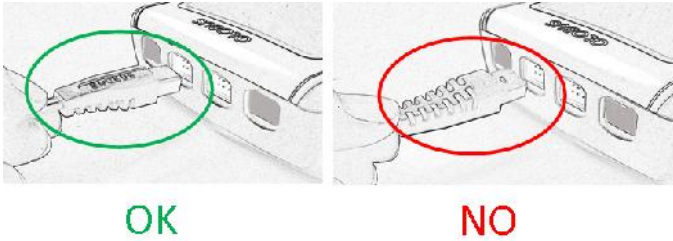E-Mat 100

The E-Mat 100 can be considered an expansion of Mat 100, since one can be connected to the other by two small linking cables. Thus the magnetic field produced by each channel of the Magnum device will be funneled on two mats and will increase the treated surface significantly.

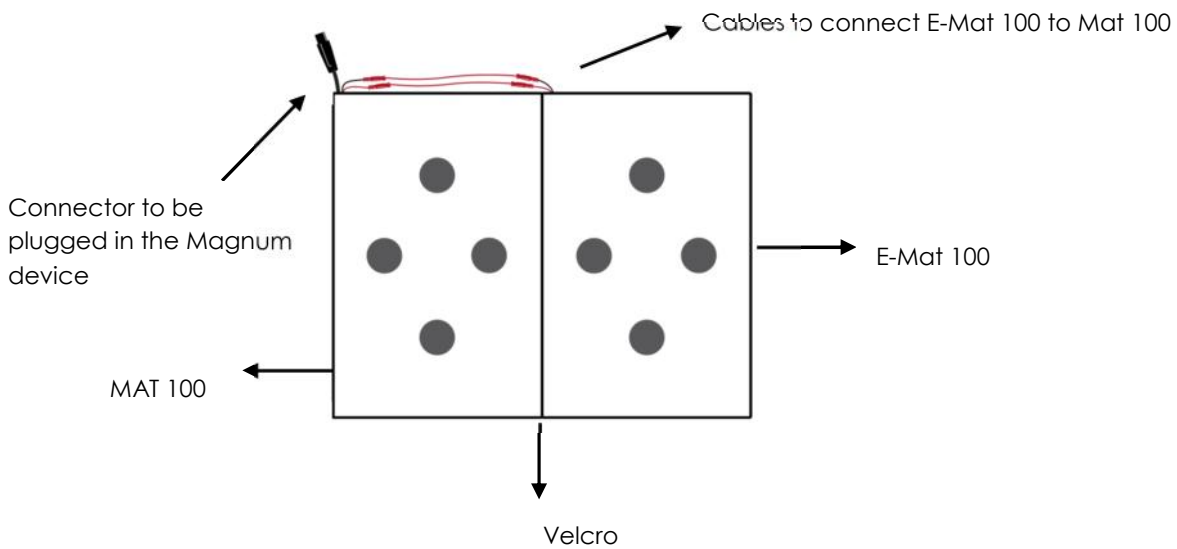To connect the two mats:


Open the Mat 100 on the shorter side and pull out the cable as shown in the picture

Unplug one of the two connectors closing the circuit

Connect the E-Mat 100 expansion by using the Velcro and pull out the cables on the shorter side

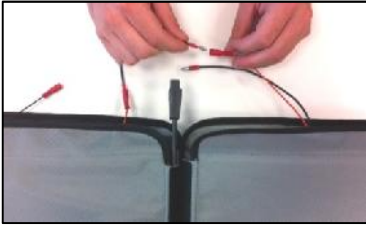Connect Mat 100 male and female connectors to the E-Mat 100 male and female ones.

Correct connection

Reintroduce all the cables inside the mat and close the Velcro

N-Mat

The N-Mat mat is a neutral mat without internal solenoids. N-Mat can be used to create a continuity support in the composition of B-Mat.

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