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Manuale

GUIDELINES AND SUMMARY ON THE DATA ACQUIRED ON THE ELECTROMAGNETIC THERAPY APPLIED BY MAGNETOMAX-PLUS

MAGNETOMAX-PLUS is a small electro-medical equipment, which allows – using its irradiating probes – to perform one of the most advanced therapies adopted until now in the leading specialized centre for the **Electromagnetic Therapy**, comfortably at home.

WHAT IS MAGNETOMAX-PLUS USED FOR?

The use of **MAGNETOMAX-PLUS** for the **Electromagnetic Therapy** is extremely simple and its harmlessness and effectiveness have been proven since a long time. Being completely harmless and painless, it can basically be used by everyone: from children and young athletes to the elderly, without having to adhere to strict limits to the duration of the therapy.

The only contraindications involve pacemaker carriers, pregnant women and people with congenital or juvenile diabetes.

Despite the very moderate price of this device (equivalent to three applications of the therapy in a specialized clinic) it meets everyone's expectations very well. In fact, it allows avoiding long waiting times at the clinic and high costs.

From this point of view, **MAGNETOMAX-PLUS** represents the most convenient and the cheapest solution.

ELECTROMAGNETIC THERAPY: A LITTLE BIT OF HISTORY

It is not long since **Electromagnetic Therapy** was officially introduced in hospitals and specialized clinics, allowing the effective treatment of the many chronic diseases that did not show signs of significant improvement with the usual pharmaceutical treatments. Despite this, the first studies on the effects of this therapy trace back to the first years of the XX Century. The first great medical school to deal with a thorough analysis on the biological effects of electromagnetic waves on body cells was the Russian one, using the first rudimentary tools.

Later, the study of some theories was successfully resumed by the French school. It was only in the Seventies that some American scholars, thanks to the new means of investigation, were able to verify that electromagnetic waves infuse the body cells with electric charge, thus improving the enzyme kinetics, speeding up tissue and bone repair, producing a healthy anti-edema, anti-inflammatory and analgesic action.

The first physicians to make the therapy official were the Germans who, treating numerous patients, were able to verify the extraordinary effectiveness of the therapy itself. They counted a conspicuous percentage of patients treated with **Electromagnetic Therapy** who obtained very satisfactory healing results. As a matter of fact, the beneficial effect of the electromagnetic waves consists in the electro-cellular recharge which rapidly re-polarizes the membranes of the damaged cells, independently of the reason (infective, chemical, mechanical or physical), restoring the normal functions.

HOW MAGNETOMAX-PLUS WORKS

Once the beneficial aspects of the **Electromagnetic Therapy** are highlighted, one might wonder how this therapy can have such a radical effect in a variety of treatments on the human body. The answer is very simple and helps to understand the principle on which this therapy is based.

The cellular membranes, especially those of the nerve cells, with their peripheral branches, have an electric potential (exactly like minuscule batteries) called “membrane at rest”.

When these cells become ill, because of an inflammatory process or other causes (among others, old age), they lose a huge part of their energy reserves, thus losing charge (again, like batteries). The discharge and the alteration of this potential are the cause of ailments which manifest in the body as inflammations, non-healing wounds, skin diseases and aching of the back, the articulations, the bone, etc.

To help all of these cells back to the normal condition, they must be recharged. Such is the goal of the **Electromagnetic Therapy**. Applied on the body part to be treated, it irradiates high-frequency electromagnetic impulses and it rapidly restores the electric potential of the discharged cellular membranes. This procedure helps and accelerates noticeably the healing process and the disappearing of localized pains.

In fact, the **MAGNETOMAX-PLUS** is equipped with two radiating solenoids that irradiates the high-frequency electromagnetic impulses.

USEFUL ADVICE ON THE USE OF THE IRRADIATING SOLENOIDS

The solenoids supplied with **MAGNETOMAX-PLUS** allows for an easy treatment. To perform the therapy, place the solenoid on the body part to be healed, then regulate the power of the channel to the desired intensity and turn the machine on for 30 - 60 minutes per session.

In the case of a cervical osteoarthritis, it is sufficient to place the probe on the neck and turn the machine on.

In the case of a fracture to the leg or the arm, it is sufficient to place the probe on the wounded part. In the case of skin diseases (acne, herpes zoster, psoriasis, eczemas, etc.) place the solenoid on the body part and perform the therapy.

The penetrating action through the body tissues of the electromagnetic waves irradiated by **MAGNETOMAX-PLUS** is around 20 cm. Thus, it is possible to perform the therapy even through bandages (in case of wounds, burns or sprains), plaster casts (in case of fractures) or clothing (in case the patient does not want to expose a specific body part). However, it is recommended to place the solenoid directly on the skin to take better advantage of the irradiating waves.

SOME THERAPEUTIC INDICATIONS

MAGNETOMAX-PLUS can treat a variety of ailments such as:

TRAUMATIC INJURIES – the treatment quickly relieves the pain and accelerates the bone consolidation. The treatment re-calcifies broken bones, cures femoral osteoporosis and fractures to the kneecap, the knee, the neck and vertebrae, thus reinforcing the bone tissue.

The therapy, in these cases, should be performed from 1 up to 3 months.

ARTHRO-RHEUMATIC DISEASES – the treatments effectively contrast muscle pain, low back pain, arthritis, osteoarthritis and a variety of bone pains (knee, hip, shoulders, elbow, heel, hands, etc.) and cervical.

After about 20 - 30 hour-long applications, the first results begin to appear and the local pain will mitigate or disappear.

CARDIOVASCULAR DISEASES – the treatments cure phlebitis, edemas, varicose veins, heavy legs, ulcers and vasculopathy.

ARTHROLYSIS AND TENDINITIS – all forms of arthrolysis and sprains (to the knee, elbow, wrist, meniscus, etc.) can be treated.

SKIN DISEASES – the treatments are very effective against acne, herpes Zoster, psoriasis, burns and eczemas.

DELAYED HEALING OF WOUNDS – the therapeutic action noticeably accelerates the healing of wounds, sores or burns. If applied on long surgical scars, after a suitable treatment, they tend to diminish. That is because the electromagnetic therapy removes dead cells and scales acting on the hemoglobin and regenerating the cells on wounds.

UROLOGICAL AND GYNECOLOGICAL DISEASES - it cures hypertrophy and prostatic inflammations, painful and irregular menstrual cycles, hemorrhoid itch.

EASE OF THE TREATMENTS

The advantages offered by the electromagnetic therapy, compared with other treatments (such as acupuncture) are undeniable. The therapy is extremely easy to perform and does not require any particular search for specific spots on the body to work on: just apply the solenoid on the body part to heal.

INTERESTING FACTS

The electromagnetic therapy has an effect on cellular regeneration and improves the general state of health, enhancing the natural defenses of the body. It has been demonstrated that its superficial action, for example, makes the skin more elastic and fresher. On a circulatory level, blood purification with a reduction of blood cholesterol has been noticed.

Good results were also achieved in the treatment of hemorrhoids with the elimination, in a short time, of itches and burns. A suitable treatment can help healing from bronchitis. What is more, the treatments have a beneficial calming effect on the nervous system.

In the USA, where the electromagnetic therapy is already widespread in gyms to fight ligament inflammations, tendinitis and to tone up muscles, many football, and basketball teams undergo intensive treatment before important games.

All medical conventions dedicated to electromagnetic therapy highlighted its extraordinary effectiveness and rapidity. Moreover, wide and extended surveys demonstrated its innocuousness and the absence of side effects.

A considerable part of the treated subjects gets a complete and satisfactory healing, the remaining part gets noticeable improvements. 10-20 60-minutes long sessions are enough to alleviate most of the localized pains. Only in case of severe inflammations or bone fractures the treatments must be prolonged up to 3 months.

Basically, positive results or improvements are always achieved.

The extreme ease of use, the absolute innocuousness, the noticeable effectiveness makes **MAGNETOMAX-PLUS** an important electric device and a valuable friend of your health **available to all**.

INSTRUCTIONS

MAGNETOMAX-PLUS provides treatment of the following pathological diseases:

- INJURIES TO THE MUSCULOSKELETAL SYSTEM AND FRACTURES
- MUSCLE PAIN, ARTHRO-RHEUMATIC DISEASES AND INFLAMMATIONS
- SKIN CONDITIONS (PSORIASIS, ACNE, ECZEMA, ALLERGIES, SORES, BURNS)
- HELPS HEALING OF WOUNDS
- INCREASES THE NATURAL DEFENSES OF THE BODY

Just a few applications are enough to alleviate pain caused by:

- ARTHRITIS, RHEUMATISMS, SPRAINS, TENDINITIS, MUSCLE PAIN
- MENSTRUAL CYCLE, PROSTATITIS, SCIATICA, NECK ARTHRITIS, MIGRAINES AND INSOMNIA.

INSTRUCTIONS OF USE

Place irradiating solenoids around the body part to treat (even if plastered or bandaged) or place directly on the shoulders, stomach, abdomen, etc. and turn the device on for at least 30-60 minutes, adjust the intensity of the channels, the device can also be used with just one solenoid. It is advisable to repeat the treatment twice a day between 8-hour intervals.

It is always best to get medical advice on the use of the magneto-therapy for every specific case, especially in case of chronic or post-surgical diseases.

CONTRAINDICATION

**Pacemaker carriers, pregnant women
and people with congenital or juvenile diabetes.**

PAUSE IMPULSE SETTING

To optimize the setting of the therapy it is important to adjust the pause impulse (knobs with flashing lights).

The numbered notches indicate the frequency of impulse:

1 - 2	For chronic diseases
3	For generic treatments
4 - 5	For acute symptoms

It is also advisable to follow this rule:

From the first to the third session, set the frequency of the impulse on high (Pos. 1-2)

From the fourth to the sixth set the frequency on position nr. 3, if the pain is not reduced, from the **seventh** session set the frequency on position 4 – 5.

The intermediate position 3 is suitable in every case without contraindications.

For prolonged therapies, such as bone fractures, treat as follows:

From the first to the third treatment set on position 2

From the fourth to the tenth set on 3

From the eleventh onwards set on position 4 – 5

The continued adjustment of the impulse, through the potentiometer, allows to lightly adjust the frequency, gradually nearing the notch recommended for the therapy. Choosing a setting over another does not cause any collateral effects but may require a longer treatment.

SPECIFICATIONS

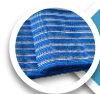
DIMENSIONS	13cm x 19cm x 6cm
WEIGHT	550 grams
TOTAL GAUSS POWER	60 Gauss (30 gauss per channel)
INTENSITY REGULATION	Set by the knobs on front
OUTPUT CHANNELS	2 Output
TREATMENT FREQUENCY	9MHz ÷ 250MHz
POWER SOURCE	220V supply
DIFFUSER	2 Solenoids
DIFFUSER DIMENSIONS	12cm x 12cm

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