mq Perfect

by maniquick®

AIR MASSAGER

**MANUALE D'ISTRUZIONI
MODE D'EMPLOI
GEBRAUCHSANLEITUNG
OPERATING INSTRUCTION
РУКОВОДСТВО ПОЛЬЗОВАТЕЛЯ**

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WARNINGS

Basic safety precautions should always be followed when using electrical appliances, especially when children are present, including the followings:

1. Keep away from water since the electrical parts are electrically live even when the power is off.
2. Always turn off immediately after use. Please remove the plug immediately if the appliance falls into water.
3. Do not place or store the device where it can fall into the water or other liquid such as the tub, sink & lavatory etc.
4. To reduce risk of burns, electrocution, fire or injury to persons, the device should never be left unattended when plugged in and in operation.
5. Under the supervision if the air leg massager is used by, on or near children or invalids.
6. Do not use accessories not recommended by the manufacturer. And should be always follow the described in this manual.
7. The device must be kept out of reach of children.

PRECAUTIONS DO NOT USE WHEN:

1. DO NOT use when bathing.
2. DO NOT use if it has a damaged cord or plug, not working properly and dropped, damaged or dropped into water.
3. DO NOT use when sleepy or not wide-awake.

This is not a medical device! This unit is to simply massage legs and feet. No medical claims are implied or warranted by the use of this product. It is strongly advised that you consult your physician before use if you are currently undertaking medical treatment.

How does it work?

Control device pumps air to gently massage legs and feet. Always have control device in hand when use. Always start with lowest power setting after that adjust the control the level what you feel comfortable.

The air leg massager requires 4xAA batteries (not included) – recommend the use of alkaline batteries for better result, or the supplied adaptor. Do not use with and extension cord only plug into the right voltage.

THE SET INCLUDING:

- A. 2 Leg boots massager
- B. 2 Thigh massager
- C. Control box
- D. 1 set of connecting pipes
- E. 6V AC/DC Adaptor


OPERATION

The air leg massager is very easy to use. You may use it either sitting or laying on a couch. With preset program to inflates and deflates automatically which can massage your legs to relieves stress and soothing fatigue. Provides a comfortable and compression massage. Gently and softly relaxes tension and improve blood circulation.

1. Outfit yourself with LEG BOOTS MASSAGER (A) and THIGH MASSAGER (B) in the applicable places. Make sure that the zipper enclosures for the boots are positioned all the way up.
2. Connect A, B and C by CONNECTING PIPES. Make sure PIPES are tightly secured to A & B.
3. Connect the AC/DC adaptor to Control box. Then plug into an outlet but make sure that the control box is OFF before plug the adaptor into an outlet.
4. Make sure to use the same type batteries (that means do not mix old with new batteries or mix alkaline with not-alkaline batteries) if you operated by batteries.
5. Press "POWER" of control box to turn on the power then the power indicator LED will bright. When you press again then the power will turn off.
6. Press "UP" or "DOWN" to adjust the working time.

IMPORTANT: ALWAYS REMEMBER THAT THE PROGRAM MUST BE ADJUSTED ACCORDING TO INDIVIDUAL SENSIBILITY. FOR THE FIRST TREATMENTS, WE ALWAYS RECOMMEND SETTING THE PROGRAM AT THE MINIMUM; THIS CAN BE GRADUALLY INCREASED ACCORDING TO YOUR SENSIBILITY, BUT ONLY THROUGH DIRECT EXPERIENCE.

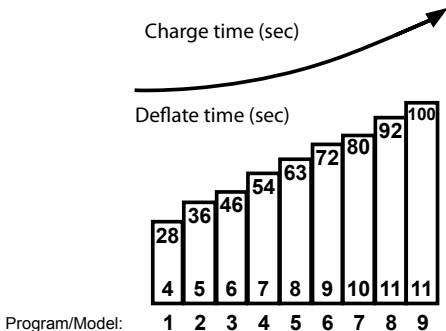IMPORTANT NOTES

Do not wash the air leg massager in wash machine. Do not dry clean. Do not bleach. Gentle wipe the device with wet cloth to gently remove the surface dust.

CONDITIONS FOR SERVICE UNDER GUARANTEE

1. SANICO Srl guarantees the absence of faults and mechanical defects and the completeness of the product when it is sold.
2. Repairs can be done only at Service Centres authorized by SANICO Srl or by the Authorized Distributor in the country where you have purchased the appliance.
3. The duration of the period of guarantee on the materials and on manufacturing and operating defects is 2 years with effect from the date of sale, subject to all the conditions for the Purchaser shown in the Instructions for use. The obligations of guarantee do not cover expendable materials: electrode cups, elastic bands, power supply cables etc..
4. Complaints relative to the product must be presented directly to the authorized service centre. The service centre, on the basis of an agreement with the Vendor, is competent to establish the cause and nature of the fault. The product will be replaced only if it cannot be repaired.
5. The owner will provide to deliver the product to the after-sales service center, as per the delivery instructions that can be requested by e-mail: info@maniquick.com.
6. The time necessary for the repair under guarantee will be subject to the degree of seriousness of the fault.
7. The responsibility of SANICO Srl under these conditions of guarantee will be limited to the amount paid by the Purchaser for the product in question.

The guarantee on the product will not be applied in the following cases

- certificate of guarantee filled in incorrectly or falsified; serial number (if applicable for the specific product) damaged or removed;
- absence of commercial and financial documents (cash receipt, delivery note) proving the purchase;
- opening or repair of the product by the Purchaser or by an unauthorized Service Centre;
- defects following failure by the Purchaser to observe the rules for use, storage and transport of the product shown in the relative technical and regulations documentation, in which case the regulations of the aforementioned documentation which have been infringed must be specified;

- penetration of foreign bodies or liquids inside the product;
- presence on the product of internal or external mechanical defects (cracks, signs of a fall, chipped parts, etc.) consequent to the incorrect use, installation or transport of the product;
- components or accessories subject to rapid wear and tear (electrode cups, elastic bands);
- defects caused by external agents (dirt visible externally and internally, traces of corrosion);
- defects due to causes independent of the manufacturer, such as power supply voltage step changes, natural phenomena or disasters, fires, actions by animals, insects etc.
- complaints for defects found by the Purchaser that are made after the expiry of the guarantee period or any other period established by the law for reporting any defects;
- use of the product by the Purchaser not for personal (domestic) use but for commercial and profit-making purposes.

Conditions for the return of a defective product

Defective products will be accepted only if clean and complete, in the original packaging, together with the certificate of guarantee correctly filled in and the commercial and financial documents (cash receipt of the store, delivery note) that prove the purchase. Please attach a description of the problem found to the appliance.

After Sales Service Center, please contact:

e-mail: info@maniquick.com

www.maniquick.com

Before sending the appliance, please contact us by phone or e-mail to agree on how the product is to be returned

SEND THIS COUPON IN CASE OF REPAIRS:

Model: MQ790

Date of purchase / /

Dealer's Stamp:

Buyer's Full Name: _____

Street/Square _____ **N°** _____

City and State: _____

Postal Code: _____

Phone Number: _____

e-mail: _____

Problem description: _____

Signature: _____

Date: _____

WARRANTY IS VALID ONLY IF ACCOMPANIED BY INVOICE/TICKET.

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DISPOSAL OF YOUR OLD APPLIANCE

When this crossed-out wheeled bin symbol is attached to a product, it means the product is covered by the European Directive 2002/96/EC.

All electrical and electronic products should be disposed off separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.

The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health.

For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

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SECURITY OF ALKALINE BATTERIES INTO PRODUCTS

1. Fire - Never delete the batteries throwing them in a fire. they could break and release the substances contained inside them.
2. If ingested or inserted, seek immediate medical assistance.
3. Do not store or leave the batteries in places accessible by children.
4. Remove old or weak batteries from appliances.
5. Do not mix old and new batteries in the same apparatus. There may be losses. The replacement of all batteries in a device must be done simultaneously.
6. Never mix different types of batteries in a single device. There may be losses.
7. Install carefully the batteries in the correct polarity (+ and -) according to the instructions provided by the manufacturer of the product
8. Do not leave the Products turn on with batteries when they are depleted. There may be losses.
9. Remove the batteries from any device that must remain unused for long periods of time
10. Avoid contact of the batteries with metal objects or other batteries because they may experience losses or they could break.
11. Do not recharge batteries because there may be some losses or they could break.
12. Not groped to disassemble the batteries
13. Store batteries in a cool and dry. Avoid sudden changes in temperature. Store the batteries in their original packaging before use.
14. RESPECT THE ENVIRONMENT - DISPOSE OF EXHAUSTED BATTERIES IN THE CONTAINERS

The battery and the possible Electronic Equipment -AEE- or anything that contains it must not be disposed as urban waste, but in compliance with Directive 2002/96/EC, 2006/66/EC and of National Laws transposition for the disposal of the product end of life, and must be subjected to a separate collection to avoid damage to the environment. Failure to comply with the above is sanctioned by Law.

IMPORTATO DA:
SANICO Srl
Via G.Ferraris, 37/41
20090 Cusago (MI) ITALY
Tel (+39) 02.90.39.00.38 - Fax (+39) 02.90.39.02.79
info@maniquick.com - www.maniquick.com